

# CD13 (Clone: EP117) Rabbit Monoclonal Antibody

PRODUCT INFORMATION:

PERFORMANCE CHARACTERISTICS:

MR1305 6ml Ready to use

MR1305 3ml Ready to use MRC1305 1ml Concentrated MRC1305 0.5ml Concentrated

MRC1305 0.1ml Concentrated MRH1305 6ml Ready to use MRH1305 3ml Ready to use

Localization: Cytoplasm / Membrane Retrieval Buffer: Tris-EDTA, pH 9.0 Incubation: 30-60 minutes

Incubation: 30-60 minutes
Positive control: Spleen, AML

#### INTENDED USE

# For research use only

This antibody is intended for use in qualitatively identify CD13 antigen by light microscopy in formalin fixed, paraffin embedded (FFPE) tissue sections using immunohistochemical (IHC) detection methodology. Interpretation of any positive or negative staining must be complemented with the evaluation of proper known controls (Positive and Negative) and must be made within the context of the patient's clinical history and other diagnostic tests. A qualified and trained pathologist must perform evaluation of the test. This antibody is intended to be used after the primary diagnosis of tumor has been made by conventional histopathology using nonimmunologic histochemical stains

#### **SUMMARY AND EXPLANATION**

CD13, also known as aminopeptidase N, was originally identified as a cell surface glycoprotein expressed by cells of granulocytic and monocytic lineages at various differentiation stages. Sequence comparisons showed that the cDNA sequence of CD13 is identical to aminopeptidase N (APN), a prominent membrane-anchored metallopeptidase expressed by the brush borders of the small intestinal and renal microvillar membrane, and also in other plasma membranes. Human APN is a receptor for one strain of human coronavirus that is an important cause of upper respiratory tract infections. Human CD13 may also mediate HCMV infection by a process that increases binding, but not its enzymatic domain. CD13 has been used as a myeloid marker. The antibody labels leukemic blasts in acute myeloid leukemia (AML) and is helpful in identifying AML subtype M0 acute lymphoid leukemia (ALL). Additionally, CD13 is a sensitive but not entirely specific marker for anaplastic lymphoma kinase positive (ALK+) anaplastic large cell lymphomas (ALCLs). CD13 is also expressed in nonhematopoietic cells including fibroblasts, bone marrow stromal cells, osteoclasts and epithelial cells. A canalicular staining pattern of CD13 in hepatocellular carcinoma (HCC) is useful in differentiation between HCC and non-HCC in liver.

# PRINCIPLE OF THE PROCEDURE

The identification of the antigen on the FFPE tissues is carried out using the above stated antibody. The antigen and antibody complex is visualized using a enzyme coupled (HRP/AP) secondary antibody with specific binding to the primary antibody, this complex is visualized by the enzymatic activation of the chromogen resulting to a visible reaction production of the antigenic site. Each and every step involves precise time and optimal temperature and the results are interpreted using a light microscope by a qualified and trained pathologist.

## REAGENT PROVIDED

Concentrated format: Antibody to CD13 is affinity purified and diluted in antibody diluent with 1% bovine serum albumin (BSA) and 0.05% of sodium azide (NaN3). Recommended dilutions: 1:50-1:100.

The antibody dilution and protocol may vary depending on the specimen preparation and specific application. Optimal conditions should be determined by individual laboratory.

Pre-diluted format: PathnSitu's ready to use antibodies are pre-tittered to optimal staining conditions. Further dilution will affect the efficacy of the antibody and may yield to sub-optimal staining.

Immunogen: A synthetic peptide corresponding to residues in human CD13 protein Host, Isotype: Rabbit, IgG

# **RUO**

#### STORAGE AND HANDLING

**Storage Recommendations:** Store at 2-8°C. When stored at the appropriate conditions, the antibody is stable until expiry. Do not use the antibody after expiration date provided on the vial in any condition.

To ensure proper regent delivery and stability, replace the dispenser cap after every use and immediately place the vial into the refrigerated conditions in an upright position. The contents of the vial should be used within 9 months from the opening of the vial.

#### SPECIMEN PREPARATION

#### Staining Recommendations:

Routinely processed, FFPE tissues are suitable for use with this primary antibody, when used PathnSitu's Poly Excel HRP/DAB detection system. The recommended tissue fixative is 10% neutral buffered formalin. Variable results may occur as are result of prolonged fixation or special processes such as decalcification. Thickness of the sections should be 2-5µm. Slides should be stained once the sections are made as antigenicity of the cut sections may diminish over a period of time. It is recommended to stain known positive and negative controls simultaneously with unknown specimens.

#### **PRECAUTIONS**

- 1. This product should be used by qualified and trained professional users only
- The product contains < 0.1% of sodium azide as preservative and is not classified hazardous, refer MSDS for further details
- As with any product derived from biological sources, proper handling procedures should be used
- Do not use reagents after expiration date
- 5. Use protective clothing and gloves, while handling reagents
- All hazardous materials should be disposed according to local state and federal regulations
- 7. Avoid microbial contamination of reagents as it may lead to incorrect results

### STAINING PROCEDURE

Antigen Retrieval Solution: Use Tris-EDTA Buffer (Cat#PS009) as antigen retrieval solution.

**Heat Retrieval Method**: Retrieve sections under steam pressure for 15 minutes using PathnSitu's MERS (Multi Epitope Retrieval System) for optimal retrieval of the epitopes, allow solution to cool at the room temperature, transfer the tissue sections/slides to the distilled water prior to the primary antibody application.

**Primary Antibody:** Cover the tissue sections with primary antibody and incubate for 30-60 min at room temperature when used PathnSitu's PolyExcel Detection System.

**Detection System:** Refer to PathnSitu's PolyExcel HRP/ DAB detection system protocol for optimal staining results.

## **QUALITY CONTROL**

The recommended positive tissue control for CD13 is Spleen and Acute Myeloid Lymphoma. A positive and negative tissue control must be run with every staining procedure performed for monitoring the correct performance of processed tissue and test reagents. A negative tissue controls provide an indication of non-specific background staining. If the results are not expected in positive and negative controls the test must be considered invalid and entire procedure must be cross verified. Individual laboratory must establish their own quality control to validate the process and antibody when opened a vial.

## INTERPRETATION OF RESULTS

CD13 stains the Cytoplasm or Membrane. A qualified experienced/trained pathologist must interpret the results in the patient's sample along with the positive and negative controls.

DS-MR1305-B. Page **1** of **2** 



**RUO** 

#### PERFORMANCE CHARACTERISTICS

PathnSitu products will undergo a thorough quality control check before it is released to the market. The antibody showed consistent specific and sensitive staining on the multiple positive tissue controls tested, by inter run, intra run and lot based studies. The antibody is stable for the expiry mentioned on the labels which is determined by real time or accelerated methods.

# TROUBLESHOOTING

- Follow the antibody specific protocol recommendations according to data sheet provided
- Tissue staining is dependent on the handling and processing of the tissue prior to staining. Improper fixation, tissue processing, antibody freezing and thawing, washing, drying, heating, sectioning or contamination with other tissues or fluids may produce artifacts, antibody trapping or inaccurate results
- 3. Do not allow the section to dry out during the entire IHC process
- Excessive or incomplete counterstaining may compromise the interpretation of the results
- If unusual results occur, contact PathnSitu's Technical Support at +91-40-2701 5544 or E-mail:<a href="mailto:techsupport@pathnsitu.com">techsupport@pathnsitu.com</a>


#### LIMITATIONS AND WARRANTY

Authorized and skilled/trained personnel only may use the product. The clinical interpretation of any test results should be evaluated within the context of the patient's medical history and other diagnostic test results. A qualified trained pathologist must perform the evaluation of the test results. There are no warranties, expressed or implied, which extend beyond the description. PathnSitu is not liable for property damage, personal injury, time or effort on economic loss caused by this product.

### **BIBLIOGRAPHY**

- 1. Bogenrieder T, et al.: Prostate 1997, 33:225-232
- 2. Shipp MA, et al.: Blood 1993, 82:1052-1070
- 3. Yeager CL, et al.: Nature 1992, 357:420-422
- 4. Soderberg C, et al.: J Virol 1993, 67:6576-6585
- 5. Scoazec JY, et al.: J Hepatol 1997, 26:543-553
- 6. Dixon J, et al.: J Clin Pathol 1994, 47:43-47
- 7. Thalhammer-Scherrer R, et al.: Am J Clin Pathol 2002, 117:380-389
- 8. Bovio IM, et al.: Am J Clin Pathol 2008, 130:628-634
- P. Rocken C, et al.: J Clin Pathol 2005, 58:1069-1075

CD13, EP117 antibody has been created by Epitomics Inc., using Epitomics proprietary rabbit monoclonal antibody technology covered under Patent No.'s 5,675,063 and 7,402,409.

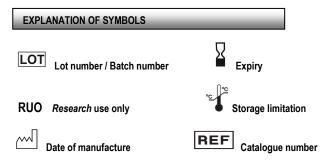

DS-MR1305-B. Page **2** of **2**